



Configuration

- Front sash with individual settings
- Light curtain for activity monitoring
- Integrated lighting
- Monitoring of filter pressure and differential pressure
- Integrated pass-through
- Monitoring of temperature
- Programmable logic controller (PLC) for monitoring and controlling function, ventilation, equipment
- Filters for supply air according ISO class

Optional configuration


- Connections for hotplates
- Bottling station (various versions)
- Underfloor flushing
- Dropping
- Media supply acc. customer requirements
- Remote maintenance
- Control units for media (hand, foot, display)
- Integrated ultrapure water system (various manufacturers)
- Still
- Ultrasonic basin
- Base cabinet with revolving doors or drawers
- Exhaust air scrubber and exhaust air irrigation

Technical data

Materials	PP, PET, POM, PFA, PVDF
Operating voltage	400 V, 5-pin
Air volume requirements	Supply air min. 400 m ³ /h, exhaust air min. 450 m ³ /h (Width: 1200 mm)
Air velocity in working area	ca. 0.15 m/s
Terminal supply air filters	High-performance HEPA filter, filter class H13 to U16 (depending on requirements) with pre-filter F9
Cleanroom classes according DIN EN ISO 14644-1	Class 3 or better (0.2 µm, 0.3 µm, 0.5 µm, 1 µm, 5 µm) „at rest“
Compliance with the DIN EN ISO 14175	Defined safety specifications of the retention and air exchange capacity for personal and product protection
Illuminance	> 750 Lux
Interfaces	Ethernet or other common fieldbus systems and serial interfaces
Dimensions	1200 x 975 x 2260 mm (W x D x H), Widths 900 mm, 1500 mm and 1800 mm also available
Weight (without equipment)	260 kg for 1200 mm width (320 kg, 380 kg for higher widths)
Dimensions working area	1180 x 690 x 880 mm (W x D x H), width varies with total width
Working height	700 mm, 900 mm
Front sash opening height	Individually adjustable via PLC



Workstation vertical

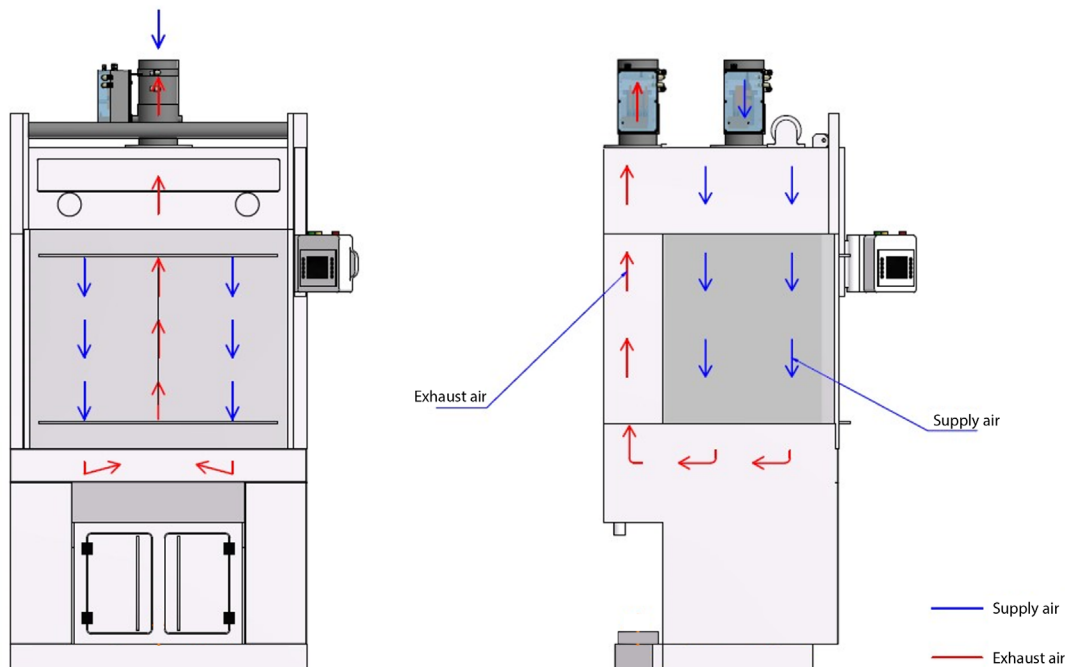

Product features

- Supply air via available supply air from in-house ventilation system
- Supply air is fed into the working area via terminal filter

- Vertical laminar flow is realized by a fine-meshed monofilament fabric

- Monofilament fabric can be changed

- Laminar air flow in the whole working area

- Control of supply air, recirculated air and exhaust air via the PLC ensures the correct pressure in the working chamber

- Compliance with DIN EN 14175

- Exhaust air can be discharged via air scrubbers or directly transferred to the exhaust air system (depending on the used chemicals)

- Room air can also be extracted via the base module and fed into the ventilation system ('recirculating air')

Contact & Support

MK Versuchsanlagen GmbH & Co. KG
Stueckweg 10 - D-35325 Muecke-Merlau
Tel.: +49 (0) 6400-9576030
info@mk-versuchsanlagen.de
www.mk-versuchsanlagen.de

